



LyteStar™ 2019-nCoV Kit

Multiplex detection of E-gene and
RdRP-gene of SARS-CoV-2


The LyteStar™ 2019-nCoV Kit is a CE IVD certified multiplex kit for detection of SARS-CoV-2 E-gene and RdRP-gene in less than 1.5 hours.

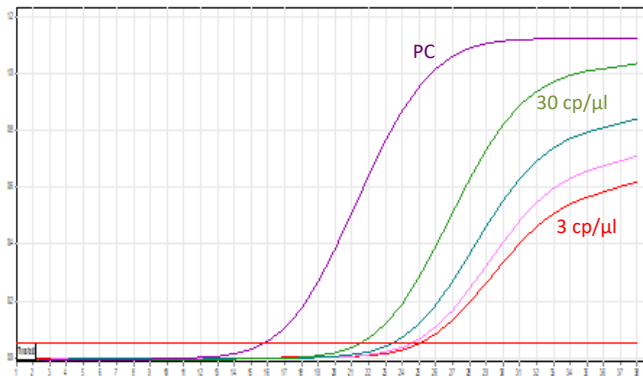


Detection on RGQ - E-gene: PC, 30, 10, 5, 3 copies/μl

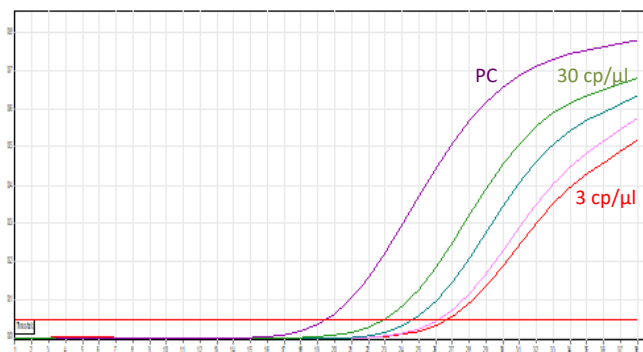

Detection on RGQ - RdRP-gene: PC, 30, 10, 5, 3 copies/μl

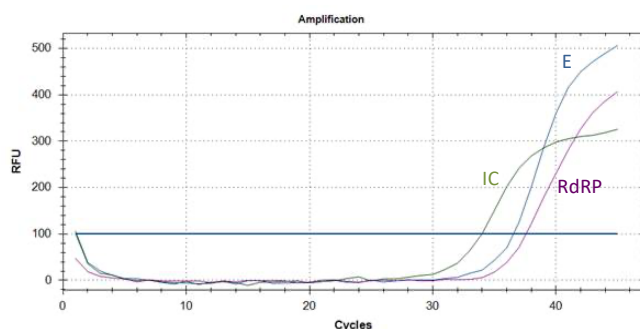

Detection on CFX96: IC, E-gene (3 cp/μl), RdRP-gene (3 cp/μl)

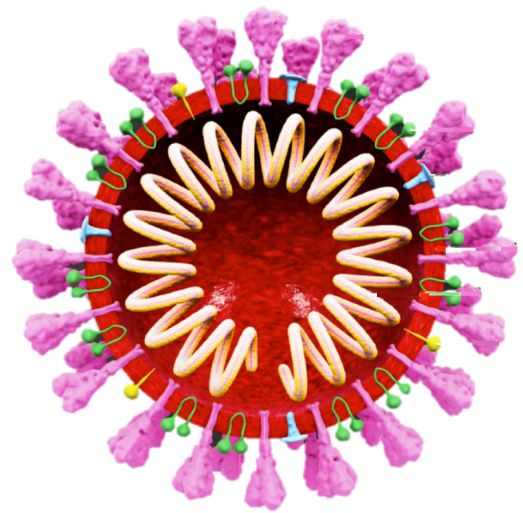

LyteStar™ 2019-nCoV RT-PCR Kit

- **Multiplex kit** for Screening and Confirmation:
 - SARS-related CoV; E-gene
 - SARS-CoV-2 specific; RdRP-gene
 - Internal Control (IC); heterologous RNA
- **Highest Analytical Sensitivity**
 - 2.14 cp/μl for E-gene
 - 1.68 cp/μl for RdRP-gene
- **100% Analytical Specificity; no cross reactivity** with other human CoV, Influenza and other respiratory viruses and pathogens
- **100% Diagnostic Sensitivity**
- **100% Diagnostic Specificity**
- For Rotor-Gene Q (Qiagen), CFX96 (BioRad), Mic (Biomolecular Systems), ABI Prism 7500 SDS and QuantStudio 5 (Applied Biosystems)
- **CE IVD certified**
- **MDA Reg. No. IVDC2905923-115842**

Recommended Sample Preparation

SpinStar™ Viral Nucleic Acid Kit

The SpinStar™ Viral Nucleic Acid Kit is designed for the simultaneous purification of viral DNA and RNA from fresh or frozen human serum and plasma, respiratory specimens and stool. The extracted DNA and/or RNA is ready to use for any downstream application where highly purified nucleic acids are required.



Validated for Direct PCR

The LyteStar™ 2019-nCoV RT-PCR Kit 2.0 has also been validated for direct PCR from Viral Transport Medium without the need for prior nucleic acid extraction. For validated starting materials refer to Instructions for Use.

Order Information

Product	Content	Σ	Cat. No.
LyteStar™ 2019-nCoV RT-PCR Kit 1.0 (<i>E</i> & <i>RdRP</i>)	Master A, Master B, Positive Control (PC), Internal Control (IC), PCR-grade Water	96	888003
SpinStar™ Viral Nucleic Acid Extraction Kit	SpinStar™ Columns, Collection Tubes, Lysis Buffer, Wash Buffer 1, Wash Buffer 2, Carrier RNA, Proteinase K, Elution Buffer	100	811803

ADT Biotech

307 Block B, Phileo Damansara 1, 9 Jalan 16/11
46350 Petaling Jaya, Selangor, Malaysia

phone +603 7931 6760
fax +603 7931 5352
email info@adt-biotech.com
website www.adt-biotech.com

